


BW-Type

Implant

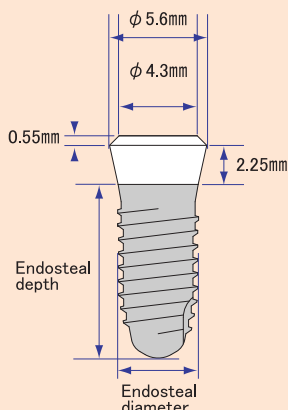

Pure Titanium
Surface treatment :
ABS


Insertion jig is
pre-mounted.

Incl. Seal Cap
[Sterilized]

The implant of this type has the same features as the B-type implant except it comes in a wider diameter (The length of a transgingival section is designed to correspond to the average thickness of the gingiva). Indications include cases requiring a back-up implantation due to the lack of initial stability after the placement of the rather shorter implant and implantation at the site of a large bone defect following extraction or loss of teeth. The BW-type implants are best suited for single tooth restorations in the posterior region of the maxilla and the mandible.

MYTIS Arrow Implant $\phi 4.7\text{mm}$

Implant (Incl. Seal Cap)

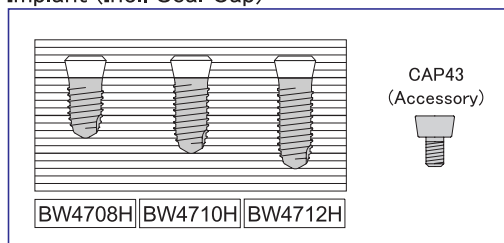

After fixture insertion, seal the internal connection of the fixture with the Seal Cap.
Installation tool for Seal Cap and Cover Screw: HEX Driver

It is used to seal the internal connection for the abutment after fixture insertion. It is used when the fixture shoulder is covered with gingiva (when the abutment is connected.)

Cover Screws

